Creative in care











Popular line of multipurpose rehab-buggies CORZO is now available as well as in Special Edition with the state of the art design outfit.

Wide range of multipurpose rehab-buggy CORZO is manufactured in 4 sizes (30, 34, 38 and 42). Buggy is determined for a short-stay usage of children with special needs, e.g. during the transportation, shopping, ..etc. Buggy line is based on well approved fully foldable undercarriage of specific authentic design, which is produced from lightweight alloy profiles with powder coating finish. Chassis is equipped with optimized buggy-folding system, wheels of 200 / 260 mm, light frame, 5-point fixation ensuring high safety and comfortable handling, high maximum load (up to 75 kg - sizes 38, 42) and very easy folding. Those features make the products very well suitable for in-town conditions and public transportation, but as well as for light outdoor activities. Despite all the advantages of the buggy-folding, CORZO Xcountry – products line preserves as well as the positioning and the "growing" advantages of the standard rehab-strollers, such a seat depth adjustment, height adjustable backrest, height adjustable footrest and backrest height extension (size 38, 42 only). Wide variety of sizes and modifications ensure fulfilling of majority requirements of user's wide spectrum. The correct usage of the product should be consulted with a seating specialist. The products corresponds to the newest EC standards as well as hygienic, technical and safety requirements.

Standard equipment:

- seat depth adjustment
- height adjustable backrest
- fully height adjustable and foldable footrest
- backrest angle adjustment (2 positions)
- backrest height extension (size 38, 42 only)
- 5-point safety harness with crotch belt and shoulder padding
- 200 mm front / 260 mm rear solid wheels
- · cover patterns and frame colours according to current offer

Optional equipment:

- hood
- comfort seat cover
- 2-piece flip-up footrest
- headrest, abductor, side trunk rest, fixation waistcoat, H-style harness, foot fixation, and other positioning accessories according to the actual list

Wide variety of fixation devices and accessories as an optional equipment (ask for accessory actual list and leaflet)

Standard foldable footrest or 2-piece flip-up footrest on request.



SUCCESSFULLY CRASH TESTED

The product had been successfully crash tested and fulfilled all the requirements of EN ISO 7176/19. It enables to use the rehab-buggy not only for a common duty, but as well as for the transportation in specially adjusted motor vehicles, equipped with the anchoring points and equipment.



Order codes - solid tyres (ST):

CORZO Xcountry-SE 30 CRX30PYYY-SE CORZO Xcountry-SE 34 CRX34PYYY-SE CORZO Xcountry-SE 38 CRX38PYYY-SE CORZO Xcountry-SE 42 CRX42PYYY-SE

Basic technical data - CORZO Xcountry-SE (CRX-SE)					
	CRX 30	CRX 34	□ CRX 38	CRX 42	
Backrest height [mm]	540-630	580-670	670-930	670-930	
Seat depth [mm]	260-340	300-370	320-390	320-400	
Seat width [mm]	300	340	1 380	420	
Seat to footrest [mm]	240-370	280-450	290-470	330-500	
Backrest angle	106°, 111°	i 110°, 114°	i 110°, 114°	i 108°, 112°	
Seat angle	17°	! 22°	21°	19°	
Total length [mm]	980	1030	1150	1150	
Total width [mm]	530	I 580	I 630	I 650	
Total height /with handle [mm]	1010-1120	980-1080	1020-1180	1210-1320	
Brake	foot oper.	foot oper.	foot oper.	foot oper.	
Weight	12,5 kg	! 13 kg	! 14 kg	! 14,7 kg	
Maximum load	45,0 kg	55,0 kg	75,0 kg	75,0 kg	
Folded stroller length [mm]	1170	i 1170	1200	1260	
Folded stroller width [mm]	320	320	350	370	
Folded stroller height [mm]	350	370	380	350	





CORZO Xcountry





Newly redesigned undercarriages in 4 sizes are the base for the CORZO's very popular line of multipurpose rehab-buggies. Wide range of rehab-buggies CORZO X country in all sizes 30, 34, 38 and 42 is determined for a short-stay usage of children with special needs, e.g. during the transportation. CORZO Xcountry - line is based on redesigned undercarriage of well approved previous version. Fully foldable undercarriage of specific authentic design is produced from lightweight alloy profiles with powder coating finish. Buggy is equipped with optimized buggy-folding system, wheels of 175 / 260 mm, light frame, 5-point fixation ensuring high safety and comfortable handling, high maximum load (up to 75 kg - sizes 38, 42) and very easy folding. Those features make the products very well suitable for in-town conditions and public transportation, but as well as for light outdoor activities. Despite all the advantages of the buggy-folding, CORZO Xcountry – products line preserves as well as the positioning and the "growing" advantages of the standard rehab-strollers, such a seat depth adjustment, height adjustable backrest, height adjustable footrest and backrest height extension (size 38, 42 only). Wide variety of sizes and modifications ensure fulfilling of majority requirements of user's wide spectrum. The correct usage of the product should be consulted with a seating specialist. The products corresponds to the newest EC standards as well as hygienic, technical and safety requirements.

Standard equipment:

- · seat depth adjustment
- height adjustable backrest
- fully height adjustable and foldable footrest
- backrest angle adjustment (2 positions)
- backrest height extension (size 38, 42 only)
- 5-point safety harness with crotch belt and shoulder padding
- 175 / 260 mm solid wheels
- cover patterns and frame colors according to current offer

Optional equipment:

- hood
- · comfort seat cover
- 2-piece flip-up footrest
- headrest, abductor, side trunk rest, fixation waistcoat, H-style harness, foot fixation, and other positioning accessories according to the actual list



wide variety of fixation devices and accessories as an optional equipment (ask for accessory actual list and leaflet)



Standard foldable footrest or 2-piece flipup footrest on request



SUCCESSFULLY CRASH TESTED

The product had been successfully crash tested and fulfilled all the requirements of EN ISO 7176/19. It enables to use the rehab-buggy not only for a common duty, but as well as for the transportation in specially adjusted motor vehicles, equipped with the anchoring points and equipment.



Order codes - solid tyres (ST):

CORZO Xcountry 30 CRX30PYYY
CORZO Xcountry 34 CRX34PYYY
CORZO Xcountry 38 CRX38PYYY
CORZO Xcountry 42 CRX42PYYY

Basic technical data - CORZO Xcountry (CRX)					
	CRX 30	CRX 34	CRX 38	□ CRX 42	
Backrest height [mm]	540-630	580-670	670-930	670-930	
Seat depth [mm]	260-340	300-370	320-390	320-400	
Seat width [mm]	300	1 340	1 380	420	
Seat to footrest [mm]	240-370	280-450	290-470	330-500	
Backrest angle	106°, 111°	ı 110°, 114°	i 110°, 114°	ı 108°, 112°	
Seat angle	17°	! 22°	! 21°	! 19°	
Total length [mm]	980	1030	1150	1150	
Total width [mm]	530	I 580	I 630	I 650	
Total height /with handle [mm]	1010-1120	980-1080	1020-1180	1210-1320	
Brake	foot oper.	foot oper.	foot oper.	foot oper.	
Weight	12,5 kg	! 13 kg	! 14 kg	! 14,7 kg	
Maximum load	45,0 kg	55,0 kg	75,0 kg	75,0 kg	
Folded stroller length [mm]	1170	i 1170	1200	ı 1260	
Folded stroller width [mm]	320	320	350	370	
Folded stroller height [mm]	350	370	380	350	







CORZINO Xcountry-SE

Comfortable line of popular rehab-buggies CORZINO Xcountry newly as an art member of Special Edition design outfit. The CORZINO's range of undercarriages are produced from lightweight profiles in 4 sizes 30, 34, 38 and 42 with the high maximum load (up to 75 kg – size 38,42) and equipped with optimized easy buggy-folding system. The new re-designed soft padded ergonomic seat cover with 5-point immobilisation belt ensure comfortable seating. Wheels of 200/260 mm, light frame, comfortable handling also make the buggy suitable for in-town conditions as well as for light outdoor activities. Wide variety of sizes and modifications, basic "growing" features – fully height adjustable footrest and detachable backrest height extension (size 38,42 only), ensure fulfilling of majority requirements of user's wide spectrum.

The correct usage of the product should be consulted with a seating specialist. The products corresponds to the newest EC standards as well as hygienic, technical and safety requirements.

Standard equipment:

- well padded ergonomical seat cover
- comfort side protection
- fully height adjustable and foldable footrest
- backrest angle adjustment (2 positions)
- 5-point safety harness with crotch belt and shoulder padding
- backrest height extension (size 38,42 only)
- 200 / 260 mm solid wheels
- · cover patterns and frame colors according to current offer

Optional equipment:

- canopy
- 2-piece flip-up footrest
- headrest, abductor, side trunk rest, fixation waistcoat, H-style harness, feet fixation, and other positioning accessories according to the actual list

wide variety of fixation devices and accessories as an optional equipment (ask for accessory actual list and leaflet)

Order codes - solid tyres (ST):

CORZINO Xcountry-SE 30 CNX30PYYY-SE
CORZINO Xcountry-SE 34 CNX34PYYY-SE
CORZINO Xcountry-SE 38 CNX38PYYY-SE
CORZINO Xcountry-SE 42 CNX42PYYY-SE



	INO Xcountry SE (CNX-SE 30	CNX-SE 34	CNX-SE 38	CNX-SE 42
Backrest height [mm]	530	670	860	860
Seat depth [mm]	270	320	350	350
Seat width [mm]	300	340	! 380	420
Seat to footrest [mm]	240-370	280-450	290-470	330-500
Backrest angle	106°, 111°	i 110°, 114°	110°, 114°	108°, 112°
Seat angle	17°	22°	21°	19°
Total length [mm]	980	1030	1150	1150
Total width [mm]	530	I 580	I 630	I 650
Total height /with handle [mm]	1010	1050	1100	1210
Brake	foot oper.	ı foot oper.	foot oper.	i foot oper.
Weight	11,5 kg	! 12,5 kg	! 13,5 kg	14,5 kg
Maximum load	45,0 kg	55,0 kg	75,0 kg	75,0 kg
Folded stroller length [mm]	1170	I 1170	ı 1200	ı 1260
Folded stroller width [mm]	320	320	350	370
Folded stroller height [mm]	350	370	380	380





CORZINO Xcountry







CORZINO Xcountry 30

CORZINO Xcountry 34

CORZINO Xcountry 38

CORZINO Xcountry 42

Comprehensive range of undercarriages in 4 sizes is the base for the CORZINO's standard line of rehab-buggies. This products line is produced in sizes 30, 34, 38 and 42 and it is determined for short-stay usage of children with special needs, e.g. during the transportation. CORZINO Xcountry product line is based on redesigned undercarriage of well approved CORZO products line. Fully foldable undercarriage of specific authentic design is produced from lightweight profiles with powder coating finished. Buggy is equipped with optimized buggy-folding system, newly redesigned well padded ergonomical seat cover, wheels of 175 / 260 mm, light frame, 5-point fixation ensuring high safety and comfortable handling, high maximum load (up to 75 kg - size 38, 42) and very easy folding. Those features make the products very well suitable for in-town conditions and public transportation, but as well as for light outdoor activities. Despite all the advantages of the buggy-folding, CORZINO Xcountry - line preserves the positioning advantages of the standard rehab-strollers. Buggy offers as well as basic "growing" features - fully height adjustable footrest and detachable backrest height extension (size 38, 42 only). Wide variety of sizes and modifications ensure fulfilling of majority requirements of user's wide spectrum. The correct usage of the product should be consulted with a seating specialist. The products corresponds to the newest EC standards as well as hygienic, technical and safety requirements.



Order codes - solid tyres (ST):	Orde	r codes	- solid to	vres ((ST):
---------------------------------	------	---------	------------	--------	-------

CNX30PYYY
CNX34PYYY
CNX38PYYY
CNX42PYYY

Standard equipment:

- well padded ergonomical seat cover
- comfort side protection
- fully height adjustable and foldable footrest
- backrest angle adjustment (2 positions)
- 5-point safety harness with crotch belt and shoulder padding
- backrest height extension (size 38 only)
- 175 / 260 mm solid wheels
- · cover patterns and frame colors according to current offer

Optional equipment:

- hood
- 2-piece flip-up footrest
- headrest, abductor, side trunk rest, fixation waistcoat, H-style harness, feet fixation, and other positioning accessories according to the actual list

Standard foldable footrest or 2-piece flip-up footrest on request



wide variety of fixation devices and accessories as an optional equipment (ask for accessory actual list and leaflet)



Basic technical data - CORZ	CNX 30 `	CNX 34	¹ CNX 38	¹ CNX 42
Backrest height [mm]	530	670	860	860
Seat depth [mm]	270	320	350	350
Seat width [mm]	300	340	380	420
Seat to footrest [mm]	240-370	280-450	290-470	330-500
Backrest angle	106°, 111°	ı 110°, 114°	ı 110°, 114°	ı 108°, 112°
Seat angle	17°	22°	21°	19°
Total length [mm]	980	1030	1150	1150
Total width [mm]	530	I 580	I 630	I 650
Total height /with handle [mm]	1010	1050	1100	1210
Brake	foot oper.	foot oper.	foot oper.	i foot oper.
Weight	11,5 kg	! 12,5 kg	! 13,5 kg	. 14,5 kg
Maximum load	45,0 kg	55,0 kg	75,0 kg	75,0 kg
Folded stroller length [mm]	1170	ı 1170	ı 1200	ı 1260
Folded stroller width [mm]	320	320	350	370
Folded stroller height [mm]	350	370	380	380





CORZINO Classic





CORZINO Classic 30

CORZINO Classic 34

CORZINO Classic 38

CORZINO Classic 42

Comprehensive-range undercarriage in 4 sizes is the base for the CORZINO's economy line of rehab-buggies. This products line is produced in sizes 30, 34, 38 and 42 and it is determined for short-stay usage of children with special needs, e.g. during the transportation. CORZINO Classic - product line is based on redesigned undercarriage of well approved CORZO products line. Fully foldable undercarriage of specific authentic design is produced from lightweight profiles with powder coating finish. Buggy is equipped with optimized buggy-folding system, wheels of 175/260 mm, light frame, 5-point fixation ensuring high safety and comfortable handling, high maximum load (up to 75 kg - size 38, 42) and very easy folding. Those features make the products very well suitable for in-town conditions and public transportation, but as well as for light outdoor activities. Despite all the advantages of the buggy-folding CORZINO Classic - line preserves as well as the positioning advantages of the standard rehab-strollers. Buggy offers as well as basic "growing" features - fully height adjustable footrest and detachable backrest height extension (size 38, 42 only). Wide variety of sizes and modifications ensure fulfilling of majority requirements of user's wide spectrum. The correct usage of the product should be consulted with a seating specialist. All products corresponds to the standard hygienic, technical and safety requirements.



fully height adjustable and foldable footrest available also in option of Flip-up Footplates

Order codes - solid tyres (ST):

CORZINO Classic 30 CNC30PYYY
CORZINO Classic 34 CNC34PYYY
CORZINO Classic 38 CNC38PYYY
CORZINO Classic 42 CNC42PYYY

Standard equipment:

- fully height adjustable and foldable footrest
- backrest angle adjustment (2 positions)
- 5-point safety harness with crotch belt and shoulder padding
- backrest height extension (size 38, 42 only)
- big wheels 175/260 mm
- cover patterns and frame colors according to current offer

Optional equipment:

- hood
- · comfort seat cover
- comfort side protection
- 2-piece flip-up footrest
- headrest, abductor, side trunk rest, fixation waistcoat,
 H-style harness, feet fixation, and other positioning accessories according to the actual list



Basic technical data - CORZINO Classic (CNC)					
	CNC 30	I CNC 34	I CNC 38	I CNC 42	
Backrest height [mm]	530	670	860	860	
Seat depth [mm]	270	320	350	350	
Seat width [mm]	300	940	380	420	
Seat to footrest [mm]	240-370	280-450	290-470	330-500	
Backrest angle	106°, 111°	ı 110°, 114°	ı 110°, 114°	ı 108°, 112°	
Seat angle	17°	22°	21°	19°	
Total length [mm]	980	1030	1150	1150	
Total width [mm]	530	I 580	I 630	I 650	
Total height /with handle [mm]	1010	1050	1100	1210	
Brake	foot oper.	foot oper.	foot oper.	foot oper.	
Weight	11,5 kg	! 12,5 kg	13,5 kg	14,5 kg	
Maximum load	45,0 kg	55,0 kg	75,0 kg	75,0 kg	
Folded stroller length [mm]	1170	ı 1170	1200	ı 1260	
Folded stroller width [mm]	320	320	350	370	
Folded stroller height [mm]	350	370	380	380	





TOM 5 Streeter-SE



Popular line of multipurpose fully foldable and adjustable rehab-strollers TOM 5 is now available as well as in Special Edition with the state of the art design outfit.

Well approved fully foldable aluminium spring-cushioned chassis is manufactured in 3 sizes. Full aluminium seat-unit is detachable, reversible and fully size and position adjustable. It includes not only highly innovative MONOSLIDE design, but as well as tilt-in-space telescope operated system and well approved flexible system for attaching of fixation devices and positioning accessories. Seat shell is covered by fully detachable soft padded textile upholstery. The correct usage of the product should be consulted with a seating specialist. The product corresponds to the newest EC standards as well as hygienic, technical and safety requirements.

Standard equipment:

- detachable and reversible seat-unit
- height adjustable backrest
- depth and width adjustment of the seat-base.
- seat-unit with fully adjustable tilt-in-space function
- fully reclining backrest with lateral supports
- · angle and height adjustable footrest
- · angle adjustable footplate
- 200 mm front solid / 295 mm rear air wheels
- · detachable safety bar with cover
- · height adjustable ergonomic handlebar
- · upholstery with reflexive safety marks
- basket, hip rest, padded frame cover
- cover patterns and frame colours according to current offer

Optional equipment:

- · hood, rain cover
- room base
- headrest, side trunk rest, 4-point immo. belt, abductor wedge and other equipment according to the actual list.

Order codes - air tyres (AT):

TOM 5 Streeter-SE MINI:

Front swivel castor-wheels

TOM 5 Streeter-SE STD:

Front swivel castor-wheels

TOM 5 Streeter-SE MAXI:

Front swivel castor-wheels

T5SWKCMYY-SE

T5SWKCYYY-SE

T5SWKCXYY-SE



SUCCESSFULLY CRASH TESTED

The product had been successfully crash tested and fulfilled all the requirements of EN ISO 7176/19. It enables to use the rehabstroller not only for a common duty, but as well as for the transportation in specially adjusted motor vehicles, equipped with the anchoring points and equipment.



Basic technical data - TOM 5 Streeter-SE (T5S-SE)					
	MINI	STD	I MAXI		
Backrest height [mm]	440-570	510-680	550-700		
Seat depth [mm]	185-270	230-330	290-400		
Seat width [mm]	220-280	250-350	270-390		
Seat to footrest [mm]	90-260	120-360	120-390		
Backrest angle	ı 0° - 175°	ı 0° - 175°	ı 0° - 175°		
Seat angle	+45/-28°	+40/-26°	+40/-26°		
Total length [mm]	800	800	800		
Total width [mm]	600	640	680		
Total height /with adj.handle [mm]	970/1130	970/1130	990/1200		
Brake	foot oper.	foot oper.	foot oper.		
Weight	16,0 kg	18,0 kg	19,5 kg		
Maximum load	35,0 kg	40,0 kg	50,0 kg		
Folded stroller length [mm]	ı 880 I	880	ı 900		
Folded stroller width [mm]	600	640	680		
Folded stroller height [mm]	400	410	420		





TOM 5 Streeter







wide flexibility, easy adjustments, user friendly manipulation with seat-unit and undercarriage

High-tech full-aluminum chassis of the TOM 5 stroller is produced by the new high-tech ALUX production method, which takes advantage of an aluminum alloy and a titan joint hardener combination. ALUX method ensures chassis low weight and high durability.

Fully foldable chassis includes among others new well designed chassis to seat interface, safe brake-system, spring-suspension, ergonomically shaped handlebar with the height adjustment.

100% ALUX full-aluminum seat is as well as produced fully by high-tech ALUX production method. Seat-unit is detachable, reversible and fully size and position adjustable. It includes not only highly innovative MONOSLIDE design, but as well as easy tilt-in-space telescope operated system and well approved flexible system for attaching of fixation devices and positioning accessories. Seat shell is covered by fully detachable soft padded textile upholstery in fashionable design. Easy handling is ensured as well as by the stroller's fully foldable frame even with attached seat.

Stroller TOM 5 has been designed in cooperation with professional designers and seating specialists and it excels not only in design but as well as in use value – it is one of the lightest strollers in the category. Additionally 260 mm (175/260 mm) solid wheels enables usage of the product not only in town but as well as for the light outdoor activities. The correct usage of the product should be consulted with a seating specialist.

The product corresponds to the newest EC standards as well as hygienic, technical and safety requirements.





wide variety of fixation devices and accessories as an optional equipment (ask for accessory actual list and leaflet)





SUCCESSFULLY CRASH TESTED

The product had been successfully crash tested and fulfilled all the requirements of EN ISO 7176/19. It enables to use the rehab-stroller not only for a common duty, but as well as for the transportation in specially adjusted motor vehicles, equipped with the anchoring points and equipment.



folded rehab-stroller





stroller with solid tyres (ST) and front swivel castor-wheels (FSCW)

stroller with solid tyres (ST) and front fixed wheels (FFW)

stroller with air tyres (AT) and front swivel castor-wheels (FSCW)

stroller with air tyres (AT) and front fixed wheels (FFW)



3 sizes of TOM 5 rehab-stroller and each in 4 versions of chassis offers unbeatable variety of 12 basic combinations for almost each driving conditions (shopping malls, city traffic, country side causeway, snowy uneven track,)

Standard equipment:

- detachable and reversible seat-unit
- height adjustable backrest
- · depth and width adjustment of the seat-base.
- seat-unit with fully adjustable tilt-in-space function
- fully reclining backrest with lateral supports
- · angle and height adjustable footrest
- angle adjustable footplate
- 260 mm (175/260 mm) solid wheels
- · detachable safety bar with cover
- · height adjustable ergonomic handlebar
- · upholstery with reflexive safety marks
- basket, hip rest, padded frame cover
- frame colors and fabric patterns according to the current offer

Optional equipment:

- · hood, rain cover
- room base
- headrest, side trunk rest, 4-point immo. belt, abductor wedge and other equipment according to the actual list.

Order codes - solid tyres (ST):

TOM 5 Streeter MINI:

Front fixed wheels T5SYKPMYY
Front swivel castor-wheels T5SWKPMYY

TOM 5 Streeter STD:

Front fixed wheels T5SYKPYYY
Front swivel castor-wheels T5SWKPYYY

TOM 5 Streeter MAXI:

Front fixed wheels T5SYKPXYY
Front swivel castor-wheels T5SWKPXYY

Order codes - air tyres (AT):

TOM 5 Streeter MINI:

Front fixed wheels T5SYK0MYY
Front swivel castor-wheels T5SWKCMYY

TOM 5 Streeter STD:

Front fixed wheels T5SYK0YYY
Front swivel castor-wheels T5SWKCYYY

TOM 5 Streeter MAXI:

Front fixed wheels T5SYK0XYY
Front swivel castor-wheels T5SWKCXYY

fabric patterns and frame colours according to the current offer



Basic technical data - TOM 5 Streeter (T5S)				
	ı MINI	· STD	ı MAXI	
Backrest height [mm]	440-570	510-680	550-700	
Seat depth [mm]	185-270	230-330	290-400	
Seat width [mm]	220-280	250-350	270-390	
Seat to footrest [mm]	90-260	120-360	120-390	
Backrest angle	ı 0° - 175°	ı 0° - 175°	ı 0° - 175°	
Seat angle	+45/-28°	+40/-26°	+40/-26°	
Total length [mm]	760 (860*)	760 (870*)	760 (870*)	
Total width [mm]	ı 590	I 630	I 670	
Total height /with adj.handle [mm]	950/1100	950/1110	970/1180	
Brake	foot oper.	foot oper.	foot oper.	
Weight	15,5 kg	17,5 kg	! 19,0 kg	
Maximum load	35,0 kg	40,0 kg	50,0 kg	
Folded stroller length [mm]	ı 840 (850*)	ı 850 (860*)	ı 860 (870*)	
Folded stroller width [mm]	590	630	670	
Folded stroller height [mm]	360	370	370	

^{*} version with front fixed wheels (FFW)





TOM 4 Xcountry





TOM 4 Xcountry MINI

TOM 4 Xcountry STD

TOM 4 Xcountry MAXI

TOM 4 Xcountry SuperMaxi42

TOM 4 Xcountry - facelift version of very popular fully adjustable rehab-strollers for handicapped children is based on improved design of the TOM-line undercarriages with 300 mm (200/300 mm) wheels, that enables usage of the products in town and as well as for outdoor activities.

Fully foldable spring-suspensioned undercarriage is produced from lightweight thin-walled tubes with powder coating finish. Undercarriage includes new chassis to seat interface, new brake-system, ergonomically shaped height adjustable handlebar. Redesigned seat-unit is detachable, reversible and fully size and position adjustable. It includes as well flexible system for attaching of fixation devices and accessories. Seat shell is covered by fully detachable soft padded textile upholstery in fashionable design. Easy handling is ensured as well by the stroller's fully foldable frame even with attached seat. All products correspond to the standard hygienic, technical and safety requirements.





Sizes MINI, STD, MAXI are produced with front fixed wheels or swivel castor-wheels, size SuperMaxi 42 only swivel castor-wheels



SUCCESSFULLY CRASH TESTED

The product had been successfully crash tested and fulfilled all the requirements of EN ISO 7176/19.

It enables to use the rehabstroller not only for a common duty, but as well as for the transportation in specially adjusted motor vehicles, equipped with the anchoring points and equipment.



FlexiPads set

Set of side width-defining pads enables very flexible width-adjustments of the seat in accordance to the present requirements of user by maximum of 50 mm. FlexiPads set consists from side pads of different width and density, 6 pieces in total













wide flexibility, easy adjustments, user friendly manipulation with seat-unit and undercarriage





Size SuperMaxi42 reinforced version, max. load up to 75kg

equipped with workboard and room base as well ensure wide flexibility of usage





4 x 4 = 16



4 sizes of TOM 4 rehab-stroller and each in 4 versions of chassis offers wide variety of 16 basic combinations for almost each driving conditions (shopping malls, city traffic, country side causeway, snowy uneven track,)

Standard equipment:

- detachable and reversible seat-unit
- height adjustable backrest and depth adjustment of the seat-base.
- width adjustable seat by FlexiPads set (optional equipment)
 seat-unit with fully adjustable tilt in space function
- fully reclining backrest
- · angle and height adjustable footrest
- 300 mm (200/300 mm) wheels
- · detachable soft side-cushions
- · detachable safety bar with cover
- · height adjustable ergonomic handlebar
- crash-tested according to EN ISO 7176/19
- frame colour and fabric patterns according to the current offer

Optional equipment:

- hood basket FLEXI-pads set
- rain cover room base headrest
- side trunk rest, 4-point immo. belt, hip rest, abductor wedge and other equipment according to the actual list.

Order codes - solid tyres (ST):

TOM 4 Xcountry MINI:

Front fixed wheels T4XYKPMYY
Front swivel castor-wheels T4XWKPMYY

TOM 4 Xcountry STD:

Front fixed wheels T4XYKPYYY
Front swivel castor-wheels T4XWKPYYY

TOM 4 Xcountry MAXI:

Front fixed wheels T4XYKPXYY
Front swivel castor-wheels T4XWKPXYY

TOM 4 Xcountry SuperMaxi 42 (reinforced version):
Front swivel castor-wheels
T4XWKPS2Y

Order codes - air tyres (AT):

TOM 4 Xcountry MINI:

Front fixed wheels T4XYK0MYY
Front swivel castor-wheels T4XWKCMYY

TOM 4 Xcountry STD:

Front fixed wheels T4XYK0YYY
Front swivel castor-wheels T4XWKCYYY

TOM 4 Xcountry MAXI:

Front fixed wheels T4XYK0XYY
Front swivel castor-wheels T4XWKCXYY

TOM 4 Xcountry SuperMaxi 42 (reinforced version):

Front swivel castor-wheels T4XWKCS2Y

fabric patterns and frame colours according to the current offer



Basic technical data - TOM 4	MINI	I STD	ı MAXI	SuperMaxi42
Backrest height [mm]	440-610	510-650	560-670	570-670
Seat depth [mm]	240-310	270-350	320-400	320-430
Seat width [mm]	260/300	1 240/280/340	1 280/320/380	1 340/380/420
Seat to footrest [mm]	140-280	200-350	240-400	240-400
Backrest angle	0° - 180°	ı 0° - 180°	ı 0° - 160°	ı 0° - 160°
Seat angle	+/-27°	+/-27°	+/-25°	+/-25°
Total length [mm]	860	860	920	1100
Total width [mm]	610	I 610	I 650	ı 700
Total height / with handle [mm]	1010/1150	1010/1150	1030/1180	1030/1180
Brake	foot oper.	foot oper.	foot oper.	foot oper.
Weight	17,5 kg	l 18,5 kg	! 21,0 kg	! 24,0 kg
Maximum load	35,0 kg	35,0 kg	50,0 kg	75,0 kg
Folded stroller length [mm]	930	ı 930	ı 960	ı 960
Folded stroller width [mm]	610	610	650	700
Folded stroller height [mm]	550	550	630	660





TOM 4 Classic





TOM 4 Classic MINI

TOM 4 Classic STD

TOM 4 Classic MAXI

TOM 4 Classic - product line of fully adjustable rehabstrollers for handicapped children is based on improved design of the TOM-line undercarriages with 300 mm (200/300 mm) solid wheels, that enables usage of the product in town and as well as for the outdoor activities.

Fully foldable undercarriage is produced from lightweight thin-walled tubes with powder coating finish. Redesigned seat-unit is detachable, reversible and fully size and position adjustable. It includes as well new flexible system for attaching of fixation devices and accessories. Seat shell is covered by

Standard equipment:

- · detachable and reversible seat-unit
- · height adjustable backrest and depth adjustment of the seat-
- width adjustable seat by FlexiPads set (optional equipment)
- seat with fully adjustable tilt in space function
- · fully reclining backrest
- · angle and height adjustable footrest
- 300 mm (200/300 mm) solid wheels
- · detachable soft side-cushions
- · detachable safety bar with cover
- · frame colour and fabric patterns according to the current

Optional equipment:

- hood
- basket
- work board
- FlexiPads set, rain cover, headrest, side trunk rest, 4-point immo.belt, hip rest, abductor wedge and other equipment according to the actual list.

Order codes - solid tyres (ST):

TOM 4 Classic MINI:

Front fixed wheels T4CYYPMYY Front swivel castor-wheels T4CWYPMYY

TOM 4 Classic STD:

T4CYYPYYY Front fixed wheels Front swivel castor-wheels T4CWYPYYY

TOM 4 Classic MAXI:

fully detachable soft padded textile upholstery in fashionable design. Easy handling is ensured as well as by the stroller's fully foldable frame even with attached seat. All the products correspond to standard hygienic, technical and safety requirements.

wide flexibility, easy adjustments, user friendly manipulation with seat-unit and undercarriage





front fixed or swivel castor-wheels option

folded stroller







fabric patterns and frame colours according to the current offer







Basic technical data - TO	M 4 Classic (T4C)
Front swivel castor-wheels	T4CWYPXYY
Front fixed wheels	T4CYYPXYY

Basic technical data - TOM 4 Classic (T4C)			
	□ MINI	· STD	MAXI
Backrest height [mm]	440-610	510-650	560-670
Seat depth [mm]	240-310	270-350	320-400
Seat width [mm]	260/300	240/280/340	280/320/380
Seat to footrest [mm]	140-280	200-350	240-400
Backrest angle	ı 0° - 180°	ı 0° - 180°	ı 0° - 160°
Seat angle	+/-27°	+/-27°	+/-25°
Total length [mm]	990	990	1070
Total width [mm]	I 620	I 620	I 660
Total height /with adj.handle [mm]	1010	1010	1060
Brake	foot oper.	foot oper.	foot oper.
Weight	! 17,0 kg	18,0 kg	20,5 kg
Maximum load	35,0 kg	35,0 kg	50,0 kg
Folded stroller length [mm]	ı 930	ı 930	ı 960
Folded stroller width [mm]	620	620	660
Folded stroller height [mm]	550	550	630







New superb design line of Jacko stroller with Art Sign of PATRON Special Edition product family as a aluminum modern stroller for children with mild postural needs where just simplicity in use, easy handling, compact size and lightweight chassis are requested. Comfort seating on non size adjustable seat is ensured by fully detachable soft padded textile upholstery in fashion style as well as by fully reclining backrest and angle and height adjustable footrest. The design of the stroller enabling folding to the very compact size as well as to reach excellent rate of the width between seat and chassis, what also makes possible for size STD to get through of entrance 600 mm. Jacko is available in two sizes and it is determined not only for city usage but due to 295 mm low profile air rear wheels also for the light outdoor activities. The correct usage of the product should be consulted with a seating specialist. The product corresponds to the newest EC standards as well as hygienic, technical and safety requirements.

Standard equipment:

- · fully reclining backrest
- · angle and height adjustable footrest
- angle adjustable footplate
- 200 mm front solid/ 295 mm rear air wheels
- · detachable safety bar with cover
- · height adjustable ergonomic handlebar
- reflexive safety marks
- frame colours and fabric patterns according to the current offer.

Optional equipment:

- · padded frame cover
- basket
- sun shade
- headrest
- 5-point immobilization belt, abductor wedge and other equipment according to the actual list.







fabric patterns and frame according to the current offer

variety of fixations devices and accessories as an optional equipment



Order codes - solid tyres (ST):

JACKO Streeter-SE STD:

Front swivel castor-wheels

J5SWKCYYY-SE

JACKO Streeter-SE MAXI:

Front swivel castor-wheels J5SWKCXYY-SE

Basic technical data - JACKO (J5S-SE)		
	□ J5S-SE STD	J5S-SE MAXI
Backrest height [mm]	630	710
Seat depth [mm]	280	350
Seat width [mm]	<u>!</u> 340	1 380
Seat to footrest [mm]	280-380	310-410
Backrest angle	0°- 175°	ı 0°- 175°
Seat angle	20°	! 20°
Total length [mm]	770	800
Total width [mm]	ı 580	ı 640
Total height /with adj.handle [mm]	950/1150	970/1200
Brake	foot oper.	foot oper.
Weight	□ 15,0 kg	I 17,0 kg
Maximum load	40,0 kg	50,0 kg
Folded stroller length [mm]	870	960
Folded stroller width [mm]	<u>!</u> 580	I 640
Folded stroller height [mm]	370	370







JACKO Streeter





JACKO Streeter STD

JACKO Streeter MAXI

Brand new stroller Jacko is based on new and lightweight aluminum chassis, accentuating the simplicity in usage and modern pure design. The stroller is intended for children with mild postural needs where just easy handling, very compact size and lightweight chassis are requested. Comfort seating on non size adjustable seat is ensured by fully detachable soft padded textile upholstery in fashion style as well as by fully reclining backrest and angle and height adjustable footrest.

The design of the stroller enabling folding to the very compact size as well as to reach excellent rate of the width between seat and chassis, what also makes possible for size STD to get through of entrance 600 mm.

Jacko is available in two sizes and it is determined not only for city usage but due to 260 mm solid wheels also for the light outdoor activities.

The correct usage of the product should be consulted with a seating specialist. The product corresponds to the newest EC standards as well as hygienic, technical and safety requirements.



- fully reclining backrest
- angle and height adjustable footrest
- angle adjustable footplate
- 260 mm (175/260 mm) solid wheels
- · detachable safety bar with cover
- · height adjustable ergonomic handlebar
- reflexive safety marks
- frame colours and fabric patterns according to the current offer.

Optional equipment:

- · padded frame cover
- basket
- sun shade
- headrest
- 5-point immobilization belt, abductor wedge and other equipment according to the actual list.



fabric patterns and frame according to the current offer







variety of fixations devices and accessories as an optional equipment

Order codes - solid tyres (ST):

JACKO Streeter STD:

Front swivel castor-wheels

J5SWKPYYY

JACKO Streeter MAXI:

Front swivel castor-wheels

J5SWKPXYY

Basic technical data - JACKO (J5S)		
	J5S STD	J5S MAXI
Backrest height [mm]	630	710
Seat depth [mm]	280	350
Seat width [mm]	340	1 380
Seat to footrest [mm]	280-380	310-410
Backrest angle	ı 0°- 175°	0°- 175°
Seat angle	20°	! 20°
Total length [mm]	770	800
Total width [mm]	ı 580	ı 640
Total height /with adj.handle [mm]	950/1150	970/1200
Brake	foot oper.	foot oper.
Weight	□ 15,0 kg	ı 17,0 kg
Maximum load	40,0 kg	50,0 kg
Folded stroller length [mm]	870	960
Folded stroller width [mm]	580	I 640
Folded stroller height [mm]	370	370







BEN 4 Plus







BEN 4 Plus MINI

BEN 4 Plus STD

BEN 4 Plus MAXI

BEN 4 Plus SuperMaxi42

Robust and time-proven durable undercarriage with ergonomically shaped handle and 300 mm (200/300 mm) wheels along with 4 seat-sizes (MINI, Standard, MAXI, SuperMaxi 42) are the base for the standard line of rehabstrollers. BEN 4 Plus - line is equipped with fully foldable undercarriage, produced of the light weight thin-walled tubes with powder coated surface finishing. Undercarriage is reinforced and optimized for heavy duty usage. Non size adjustable seat with option of 3 angle positions and height adjustable footrest is covered by fully detachable soft padded textile upholstery in fashionable design. The stroller can be optionally equipped with wider variety of special fixing devices and accessories. Easy handling is ensured by the stroller's fully and very compact foldable frame as well by the stroller's big wheels. All products correspond to the standard hygienic, technical and safety requirements.

Standard equipment:

- non size-adjustable seat with fully detachable soft padded textile upholstery
- fully reclining backrest
- seat with tilt-in space function (3 positions)
- height variations for attaching of the fixation devices and accessories
- angle and height adjustable footrest
- 300 mm (200/300 mm) wheels
- detachable safety bar with cover
- frame color and fabric patterns according to the current offer

Optional equipment:

- · hood, rain cover
- headrest, 4-point immo.belt, abductor wedge and other equipment according to the actual list
- work board
- basket







front fixed or swivel castor-wheels option



Order codes - solid tyres (ST):

BEN 4 Plus MINI:

Front fixed wheels B4PYYPMYY
Front swivel castor-wheels B4PWYPMYY

BEN 4 Plus STD:

Front fixed wheels B4PYYPYYY
Front swivel castor-wheels B4PWYPYYY

BEN 4 Plus MAXI:

Front fixed wheels B4PYYPXYY
Front swivel castor-wheels B4PWYPXYY

BEN 4 Plus SuperMaxi 42:

Front swivel castor-wheels B4PWYPS2Y

	ı MINI	I STD	I MAXI	□ SuperMaxi42
Backrest height [mm]	530	630	710	710
Seat depth [mm]	250	280	360	380
Seat width [mm]	<u>'</u> 300	1 340	970	420
Seat to footrest [mm]	130-260	190-330	220-380	220-380
Backrest angle	ı 0° - 180°	ı 0° - 180°	ı 0° - 130°	ı 0° - 130°
Seat angle	17°/ 20°/ 23°	17°/ 20°/ 23°	17°/ 20°/ 23°	17°/ 20°/ 23
Total length [mm]	1020	1020	1070	1100
Total width [mm]	I 620	ı 620	I 660	I 700
Total height /with adj.handle [mm]	1045	1045	1060	1060
Brake	foot oper.	foot oper.	foot oper.	foot oper.
Weight	! 15,0 kg	! 16,5 kg	! 17,5 kg	! 21,0 kg
Maximum load	35,0 kg	35,0 kg	50,0 kg	75,0 kg
Folded stroller length [mm]	i 920	ı 920	930	960
Folded stroller width [mm]	610	610	650	700
Folded stroller height [mm]	580	580	530	530





Optional accessories for rehab-strollers



headrest adjustable



headwings



headrest with neckrest



lumbar support



adjustable



five-point immobilisation belt



immobilisation belt



two-point immobilisation belt



chest-belt addition to 4 & 5-point belt



immobilisation



H-style harness



H-style harness TBX (breathable material)



chest immobilisation



Tilt lock for spastic children



abduction block



abduction belt



Pelvis fixation



knee block



Calf rest with fixation



feet fixation





summer footmuff



winter footmuff





with side protection









heavy duty basket













Optional acessories for rehab-buggies









































TAMPA





TAMPA - position and size-adjustable seating system for children with special needs who require support and stabilization in sitting for everyday activities performed at school, at home or in a rehabilitation clinic. TAMPA offering many points of adjustable support that allow the chair to be adapted to each child's individual needs, using wide accessory range. The base of the seating system provides stability and mobilty with tilt-n-space adjustement and the seat with height adjustable armrest. TAMPA is available in two sizes – STANDARD and MAXI.

Standard equipment:

- · seat depth adjustment
- · height adjustable backrest
- · width adjustable seat by FlexiPads set
- · seat height and tilt in space function
- angle and height adjustable footrest
- · height adjustable armrest
- · rear castor wheels with brake
- telescopic wheeled base with length adjustment
- comes with safety bar, armrest, push handle, FlexiPads set and 1-piece headrest

Optional equipment:

workboard

Order codes:

• side trunk rest, 4-point immo.belt, abductor block and much more of other equipment according to the actual list.

Order codes.	
TAMPA STD	TMP101YYY
TAMPA MAXI	TMP101XYY





Fabric patterns and frame colours according to the current offer.

Basic technical data - TAMPA (TMP)		
	STD	MAXI
Backrest height [mm]	I 520 - 650	I 570 - 650
Seat depth [mm]	280 - 340	310 - 380
Seat width [mm]	i 240/280/340	i 280/320/380
Seat to floor [mm]	<u>!</u> 530 - 580	530 - 580
Seat to footplate [mm]	260 - 400	260 - 400
Backrest angle [mm]	ı 0° - 180°	ı 0° - 160°
Seat angle [mm]	. +/-27°	+/-25°
Total length [mm]	880	880
Total width [mm]	I 460	I 460
Total height / with handle [mm]	970 - 1020	970 - 1020
Brake	i foot oper.	foot oper.
Total weight [kg]	<u>!</u> 18	20
Maximum load [kg]	35	50





TOM 5 Streeter-NEWBORN





Tom 5 Streeter–NEWBORN - New product solution for handicapped newborns and infants combining advantages of comfortable baby carrycot AMIGO and rehab chassis of Tom 5 as high-tech full-aluminum and durable chassis which can be optionally used also with Heavy duty steel basket for respiratory equipment or even a baby car seat MIMMO Plus (gr.0+) or other equipment.

After the baby will grow up the carrycot AMIGO can be easily replaced by rehab seat unit and further used as a multipurpose fully foldable and adjustable rehab stroller. The correct usage of the product should be consulted with neonatologists and seating specialist. The product corresponds to the newest EC standards as well as hygienic, technical and safety requirements.

Standard equipment:

- baby carrycot AMIGO
- detachable CSI-interface with clip-A-go adapters and fully adjustable tilt-in-space function
- 200 mm front solid / 295 mm rear air wheels
- · height adjustable ergonomic handlebar
- basket, padded frame cover
- · frame colour and fabric pattern according to the current offer

Optional equipment:

- Car seat MIMMO Plus (gr. 0+)
- HD -basket (steel or textile)
- hand operated drum brakes
- rehab seat unit Tom 5 with full additional accessories options

Order codes - air tyres (AT): TOM 5 Streeter-NEWBORN:

Front swivel castor-wheels

T5SWKCNBY







Basic technical data - TOM 5 Streeter-NEWBORN (T5S-NB)	ı
Carrycot lenght [mm]	I 820
Carrycot width [mm]	340
Carrycot height [mm]	200
Total length [mm]	760
Total width [mm]	670
Total height /with adj.handle [mm]	ı 970/1180
Brake	foot oper.
Weight	16 kg
Maximum load with carrycot / with rehab seat unit	I 20,0/50,0 kg
Folded chassis length [mm]	860
Folded chassis width [mm]	i 670
Folded chassis height [mm]	370

Manufactured by:



DURANGO & DENVER





DURANGO and DENVER product line of telescopic chassis for outdoor (DURANGO) and indoor (DENVER) activity. Telescopic chassis are on purpose designed as a carrier for either standardised fully position and size adjustable seating system for children with special needs - TOM 5 (available in 3 sizes-MINI, STANDARD, MAXI), or as a carrier of made to measure seat-shells for seriously handicapped children with special needs (CSI-interface accessory required, not included in standard equipment). Telescopic frames are equipped with foot operated telescopic high-low adjustment and foot operated telescopic tilt in space adjustment to support and stabilize the sitting by everyday performed activities at school, at home, at rehabilitation clinic, etc. Both telescopic chassis are adjustable from 25 kg to 50 kg of user's weight, when seat unit is 5 kg. DURANGO may be optionally equipped with hand operated drum brakes. The correct usage of this product line should be consulted with a seating specialist. Products corresponds to the newest EC standards as well as hygienic, technical and safety requirements.

Standard equipment:

- · telescopic high-low adjustment function
- telescopic tilt in space adjustment function
- front and rear castor wheels with brakes (DENVER)
- 200 mm front solid / 300 mm rear air wheels (DURANGO)
- height adjustable and detachable ergonomic handlebar (DURANGO)
- frame colours (cover patterns) according to current offer

Optional equipment:

- CSI-interface, angle and height adjustable footrest set
- hand operated drum brakes (DURANGO)



high-low adjustment range



DURANGO with hand operated drum brakes

Order codes DURANGO/DENVER:			
DURANGO	DRKC161YYY		
DENVER	RKC151YYY		

Basic technical data - DURANGO/DENVER			
	DURANGO	DENVER	
Interface to floor [mm]	280-650	250-560	
Interface angle	+35/-10°	+35/-10°	
Total length [mm]	950	670	
Total width [mm]	645	510	
Total height with handle [mm]	i 840/1130	i -	
Folded frame height [mm]	475	320	
Brake	foot oper.	foot oper.	
Weight	16,0 kg	9,0 kg	
Maximum load - user + 5 kg seat unit	25,0 - 50,0 kg	25,0 - 50,0 kg	







TOM 5 Streeter - CSI MINI

TOM 5 Streeter - CSI STD

TOM 5 Streeter - CSI MAXI

Purpose made line of specially adjusted products, coming out from high-tech full-aluminum chassis of the Tom 5 stroller's line, includes the CSI-interface. This products line is strongly determine for the seat-shell producers or for the seating professionals to provide customised solution of the seating arrangements. TOM 5 CSI-version - is produced by ALUX production method, ensures chassis very low weight and high durability. Fully foldable chassis to the extra compact size, includes among others chassis to seat-schell interface with very friendly tilt in space adjustment, safe brake system, spring-suspension, ergonomically shaped handlebar with the height adjustment. Additionally 260 mm (175/260 mm) solid wheels enables usage of the product not only in town but as well as for the light outdoor activities. Products correspond to the newest EC as well as standard hygienic, technical and safety requirements.

Standard equipment:

- detachable and reversible CSI-interface with the quick-release lock
- CSI-interface with fully adjustable tilt in space function
- 260mm (175/260 mm) solid wheels
- · height adjustable ergonomic handlebar
- basket, padded frame cover
- chassis in CRASH TESTED version according to EN ISO 7176/19
- frame colour according to the current offer

Optional equipment:

- angle and height adjustable footrest set
- air wheels (rear)



folded TOM 5 STREETER CSI Version



Angle and height adjustable footrest set (optional equipment)

3 sizes of purpose-adjusted aluminium chassis with the CSI-interface

Order codes - solid tyres (ST):

TOM 5 Streeter CSI MINI:

Front swivel castor-wheels (FSCW) T5SWKPMY9

TOM 5 Streeter CSI STD:

Front swivel castor-wheels (FSCW) T5SWKPYY9

TOM 5 Streerer CSI MAXI:

Front swivel castor-wheels (FSCW) T5SWKPXY9

Chassis produced in CRASH TESTED version according to EN ISO 7176/19.

Basic technical data - TOM 5 Streeter - CSI (T5S-CSI)			
	MINI	' STD	I MAXI
Width-clearance [mm]	360	400	440
Interface to floor [mm]	i 485	ı 485	485
Interface angle	! +45°/-28°	+40/-26°	+40/-26°
Total length [mm]	760	760	760
Total width [mm]	ı 590	I 630	ı 670
Total height /with adj.handle [mm]	950/1110	950/1110	970/1180
Brake	foot oper.	foot oper.	foot oper.
Weight	11,0 kg	11,5 kg	12,5 kg
Maximum load	35,0 kg	40,0 kg	50,0 kg
Folded stroller length [mm]	i 840	i 850	i 860
Folded stroller width [mm]	590	630	670
Folded stroller height [mm]	360	370	370








Purpose-made line of specially adjusted products, which comes out from approved conception of PATRON well-known TOM-strollers' line, includes purpose-adjusted chassis and the CSI-interface. This products line is strongly determine for the seat-shell producers or for the seating professionals to provide customised solution of the seating arrangements.

TOM 4 Xcountry CSI-version - product line of fully adjustable special undercarriages is based on well approved design of the TOM-line undercarriages with 300 mm (200/300 mm) solid wheels, that enables usage of the products in town, as well as in light terrain. Fully foldable spring-cushioned undercarriage is produced from lightweight thin-walled tubes with (Footrest set available powder coating finish. Undercarriage includes on request) chassis to seat interface with the quick-release lock, new brake-system, ergonomically shaped height adjustable handlebar. Easy handling is ensured as well by the stroller's fully foldable frame. All products correspond to the standard hygienic, technical and safety requirements.

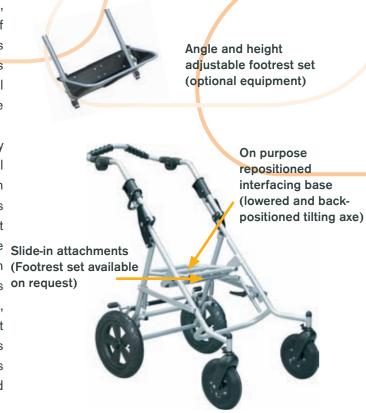

Standard equipment:

- detachable and reversible CSI-interface with the quickrelease lock
- CSI-interface with fully adjustable tilt in space function
- front fixed wheels(FFW) or front swivel castorwheels(FSCW) version
- 300mm (200/300 mm) air wheels
- height adjustable ergonomic handlebar
- chassis in CRASH TESTED version according to EN ISO 7176/19
- frame colour according to the current offer

Optional equipment:

- angle and height adjustable footrest
- air wheels (rear)

3 sizes of purpose-adjusted chassis with the CSI-interface

Order codes - solid tyres (ST):

TOM 4 Xcountry CSI-LB Standard:

Front fixed wheels (FFW) T4XYKPYT8
Front swivel castor-wheels (FSCW) T4XWKPYT8

TOM 4 Xcountry CSI-LB MAXI:

Front fixed wheels (FFW) T4XYKPXT8
Front swivel castor-wheels (FSCW) T4XWKPXT8

TOM 4 Xcountry CSI-LB SuperMaxi 42:

Front swivel castor-wheels (FSCW) T4XWKPS28

Chassis produced in CRASH TESTED version according to EN ISO 7176/19.

Basic technical data - TOM 4 Xcountry - CSI-LB (T4X-CSI-LB)			
	STD	I MAXI	SuperMaxi42
Width-clearance [mm]	370	410	460
Interface to floor [mm]	i 490	i 500	500
Interface angle	! +35°/0°	+35°/0°	+35°/0°
Total length [mm]	860	920	1100
Total width [mm]	ı 610	ı 650	700
Total height /with adj.handle [mm]	1010/1150	1030/1180	1030/1180
Brake	foot oper.	foot oper.	foot oper.
Weight	10,0 kg	11,0 kg	12,0 kg
Maximum load	40,0 kg	55,0 kg	80,0 kg
Folded stroller length [mm]	i 930	i 960	i 960
Folded stroller width [mm]	610	650	700
Folded stroller height [mm]	430	450	450





TOM 4 Xcountry DUO

TOM 4 Xcountry DUO is based on the very well approved TOM-line design of rehab-strollers for handicapped children. The fully foldable spring-suspensioned undercarriage is produced from light weight thin-walled tubes with powder coated surface finishing. Seats of the stroller are detachable and can be independently placed in or against the driving direction. Seat shells are covered by fully detachable soft padded fabric upholstery in fashionable design. Easy handling is ensured by the front swivel castor-wheels with direction lock and as well by the stroller's fully foldable frame even with attached seats. The product corresponds to standard hygienic, technical and safety requirements.

Basic equipment:

- two in function independent, detachable and independently reversible seats
- seats with fully adjustable tilt in space function
- fully reclining backrests
- · height adjustable backrest
- adjustable depth of the seats with corresponding adjustment of footrests position
- width adjustable seats by FlexiPads set (optional equipment)
- angle and height adjustable footrests
- 200mm front swivel castor-wheels / 300mm rear air wheels
- · detachable side soft cushions
- · detachable safety bars
- frame colour and fabric patterns according to the current offer

Optional equipment:

- hood
- FlexiPads set, rain cover, headrest, side trunk rest, 4-point immo.
 belt, hip rest, abductor block and other equipment according to the actual list





Order codes:

TOM 4 Xcountry DUO - MINI/MINI:

Front swivel castor-wheelsT4XWYCMDM

TOM 4 Xcountry DUO - MINI/STD:

Front swivel castor-wheels......T4XWYCCDM

TOM 4 Xcountry DUO - STD/STD:

Front swivel castor-wheels T4XWYCYDY

Basic technical data - TOM 4 Xcountry DUO			
	T4XWD-M/M	T4XWD-M/S	T4XWD-S/S
Backrest height [mm]	440-610	combi	510-650
Seat depth [mm]	240-300	ı combi	1 240-300
Seat width [mm]	240/280	combi	240/280/340
Seat to footrest [mm]	150-300	combi	200-350
Backrest angle	! 0°-180°	0°-180°	! 0°-180°
Seat angle	+/-27°	+/-27°	+/-27°
Total length [mm]	1020	1020	ı 1020
Total width [mm]	980	980	980
Total height /with adj.handle [mm]	1010	1010	1010
Weight	ı 28	ı 28	ı 28
Maximum load	60	60	60
Folded stroller length [mm]	900	900	900
Folded stroller width [mm]	ı 980	ı 980	ı 980
Folded stroller height [mm]	530	530	530
Brake	foot oper.	foot oper.	foot oper.

Manufactured by:

PATRON Bohemia a.s. | Ceskolipska 3419, CZ-27601 Melnik, Czech Republic tel.: +420-315-630-135, +420-315-630-111 | fax: +420-315-630-199 e-mail: reha.en@patron.eu | http://www.patron.eu



TOM 4 Xcountry DUO - COMBI



Tom 4 DUO-COMBI is unique product solution for siblings, offering an advanced combination of rehab seat unit with the AMIGO carrycot for newborn baby or optionally also with the ERGO sporty seat unit later on. Tom 4 DUO-COMBI is based on quality approved TOM 4 – line, fully foldable, spring-suspensioned undercarriage, equipped with rear air wheels (300 mm) and front swivel castor wheels (200 mm) ensuring easy handling. Comfortable seat unit and carrycot in fashionable design are detachable and can be placed in or against the driving direction. In addition both units can be used independently with their single original chassis to keep maximum components variability. The product corresponds to EC standard hygienic, technical and safety requirements.

Rehab seat

Basic Equipment:

- · detachable and reversible seat
- seat with fully adjustable tilt in space function
- fully reclining backrest
- · height adjustable backrest and depth adjustment of the seat-base
- width adjustable seat by FlexiPads set (optional equipment)
- · angle and height adjustable footrest
- · detachable side soft cushions
- · detachable safety bar
- frame colour and fabric patterns according to the current offer

Optional equipment:

Hood, FlexiPads set, headrest, side trunkrest, 4-point immo belt, hip rest, raincover, abductor block a other equipment according to the actual list.

Carrycot for newborn

Basic Equipment:

AMIGO carrycot with oversized width (340 mm) and length (800mm) removable upholstery with pockets and hood with summer mesh reflexive safety marks

Optional Equipment

ERGO sporty seat unit as follow-up solution summer accessory set, winter accessory set MIMMO Plus car seat

Order codes:

TOM 4 Xcountry DUO-COMBI - MINI:

Front swivel castor-wheels T4XWYCMDC

TOM 4 Xcountry DUO-COMBI - STD:

Front swivel castor-wheels T4XWYCYDC







AMIGO carrycot



MIMMO car seat

Basic technical data - TOM 4 Xcountry DUO-COMBI				
	T4XWDC-M	T4XWDC-S		
Backrest height [mm]	440-610	510-650		
Seat depth [mm]	240-300	240-300		
Seat width [mm]	240/280	240/280/340		
Seat to footrest [mm]	150-300	200-350		
Backrest angle	0°-180°	□ 0°-180°		
Seat angle	+/-27°	+/-27°		
Total length [mm]	1020	1020		
Total width [mm]	1040	1040		
Total height /with adj.handle [mm]	1010	1010		
Weight	1 24	ı 24		
Maximum load	60	<u>'</u> 60		
Folded stroller length [mm]	900	900		
Folded stroller width [mm]	1040	1040		
Folded stroller height [mm]	530	530		
Brake	foot oper.	foot oper.		



PATRON Bohemia a.s. | Ceskolipska 3419, CZ-27601 Melnik, Czech Republic tel.: +420-315-630-135, +420-315-630-111 | fax: +420-315-630-199 e-mail: reha.en@patron.eu | http://www.patron.eu



PATRON Bohemia a.s.

Ceskolipska 3419, CZ-27601 Melnik, Czech Republic

Tel.: +420 315 630 135, Fax: +420 315 630 199 E-mail: reha.en@patron.eu

http://www.patron.eu

